



GUARDIAN

SPORT

REHABILITATOR™

*The Future Standard of Care in
the World of Knee Bracing*

**IMPROVE PATIENT
OUTCOMES!**

- Increases quadriceps and hamstring strength
- Accelerates recovery
- Provides superior support, comfort and stability
- Reduces pain and inflammation
- Improves knee extension
- TREATS CHRONIC KNEE PAIN AND WEAKNESS!


1-877-371-1102 or 905-639-7243



BRACE SELECTION	INDICATION	UNLOADS	STABILITY	PAIN RELIEF	UNBRACED PAIN RELIEF AFTER 90 DAYS	IMPROVE QUADS / HAMSTRING STRENGTH	IMPROVES KNEE EXTENSION	IMPROVES FOOT PLACEMENT	GAIT RETRAINING WITH RETAINED EFFECT
SPORT REHABILITATOR™	Accelerate Recovery from Injury/Surgery	✓	✓	✓	✓	✓	✓	✓	✓
COMPETITIVE FUNCTIONAL SUPPORT BRACE	Knee Support		✓	✓					

CE U.S. Patents
7,608,051 / 7,963,933 / 8,057,414 / 8,308,669 / 8,376,947

The *SPORT* REHABILITATOR™ brace was developed for people with knee injuries including acute, chronic and post-operative conditions, as well as various instabilities. This brace offers superior support and stability to the knee as well as the ability to re-establish a normal gait pattern enhancing the rehabilitation process of the knee through a more functional neuromuscular re-education approach.

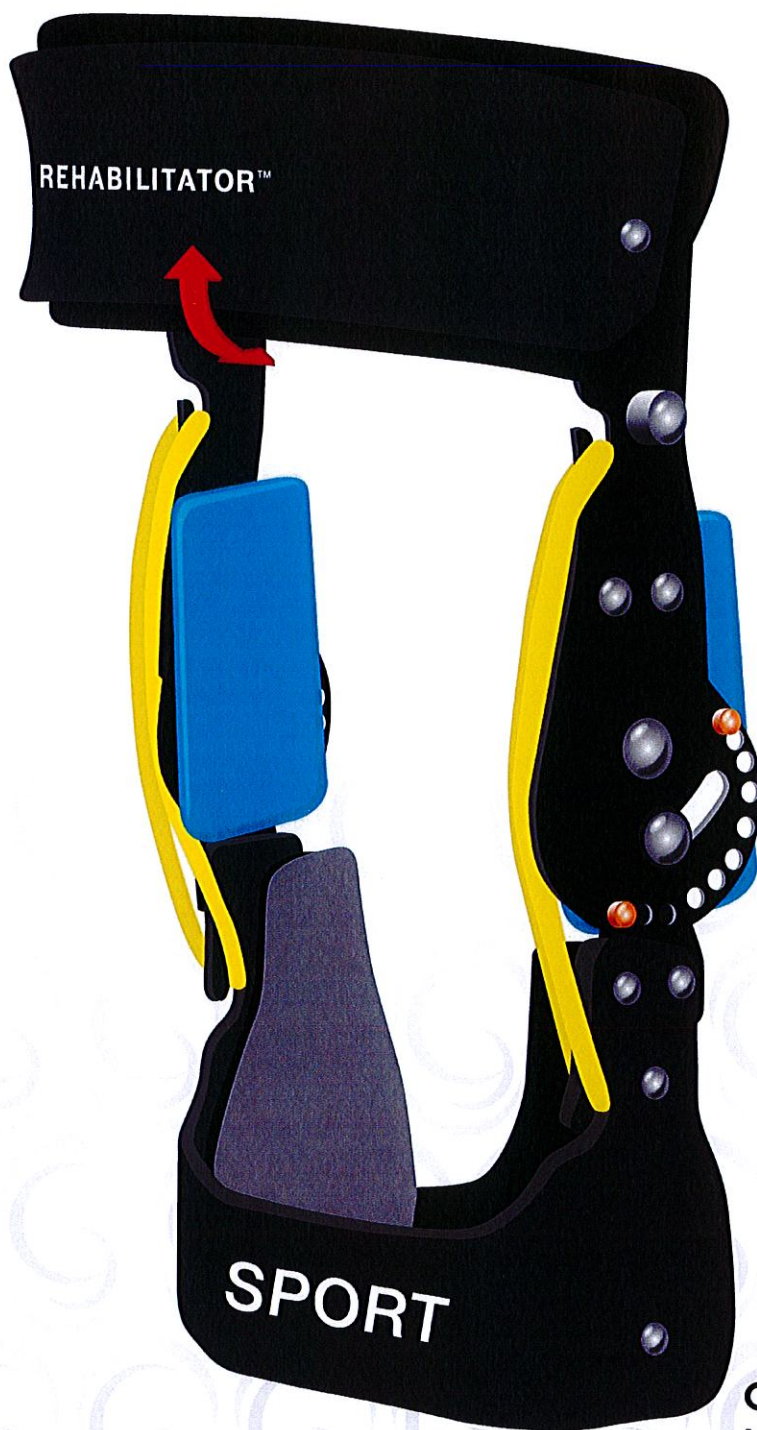

- Patented dynamic dual **SWING ASSIST™** facilitates initial swing phase of gait assisting the quadriceps with knee extension.
 - For the ACL patient with a quadriceps avoidance gait pattern, the dual **SWING ASSIST™** provides assistance in the critical open chain component of knee extension where undesired anterior translation of the tibia occurs. This can prevent or correct quadriceps avoidance gait, maintaining proper firing patterns and accelerating the post-injury or post-surgical rehabilitation process.
 - For other knee injuries that have an effusion, the dual **SWING ASSIST™** provides the same assistance with knee extension, but in this scenario the re-establishing of a normal gait pattern enhances the rehabilitation process of the quad through a more functional neuromuscular re-education approach. By assisting the quad in the open chain component of gait, the stride is lengthened, allowing a proper closed chain (eccentric loading) component and re-establishing the proper muscle firing sequence.
- Unique adjustable **pneumatic bladder system** distracts and supports the knee joint, reducing pain and inflammation while providing exceptional support.
- **ROM control hinge** allows for precise control of extension and flexion if needed.
- **Rotating thigh cuff** provides superior comfort and support.

Call 1-877-371-1102 or 905-639-7243
to order or for more information.



GUARDIAN
• BRACE •



Exclusive Canadian Distributor:
joint solutions
ALLIANCE CORPORATION
Serving Canada's Musculoskeletal Professionals

975 Fraser Drive - Unit 18 - Burlington Ontario L7L 4X8 Canada